



PocketSwab® Procedural Flow Chart using novaLUM II X


Refer to Current Charm® Operator's Manual for Complete Test Procedure

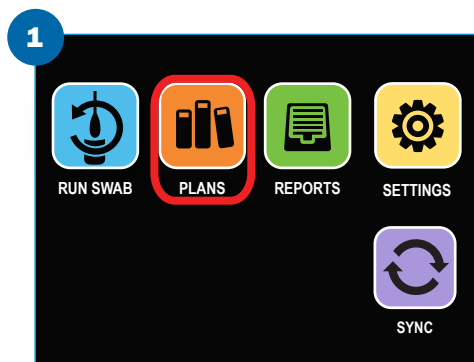


Select PLANS button.

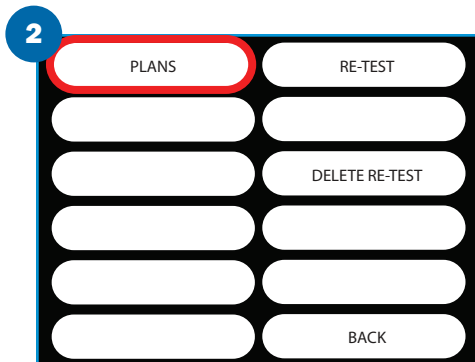

Select PLANS.

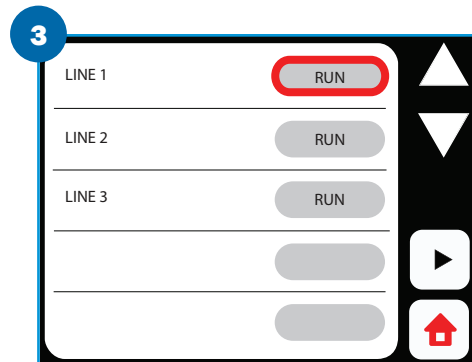

Select RUN.

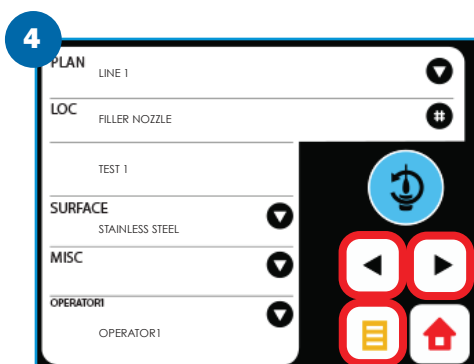

Use Left (◀) or Right (▶) arrows to pick location or use List View (☰) to scroll to find correct location.




Withdraw the swab by gently pulling and twisting the swab handle.




Swab an area 100 cm² (16 in²). Rotate the handle to achieve full swab contact with surface.



Reinsert the swab and twist the handle down. Gently shake side to side 3 times. White tablet will break apart.

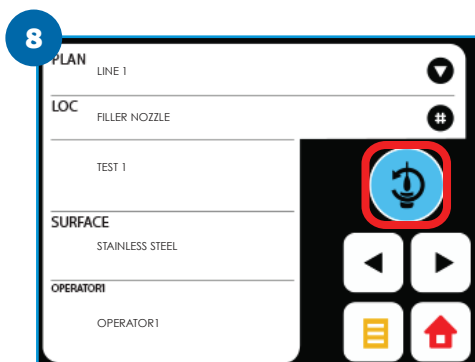

Insert the PocketSwab Plus swab fully into the chamber of the novaLUM® II system in an upright position. Press  to start count.

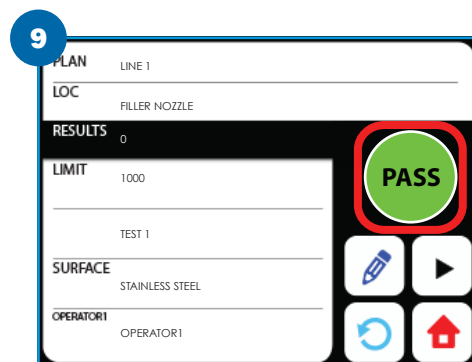

Results in 5 seconds. The result of "Pass" or "Fail" is based on established limit.

Note: Swabs can be held for up to 6 hours at room temperature before activation.

**PocketSwab®
Plus**



NOTHING WORKS LIKE A CHARM

659 Andover Street, Lawrence, MA 01843-1032, USA

T +1.978.687.9200 | F +1.978.687.9216 | E info@charm.com | www.charm.com

© 2018 Charm Sciences, Inc. Charm, novaLUM and PocketSwab are registered trademarks of Charm Sciences, Inc.

MRK-1234-002 May-2018