




WaterGiene Procedural Flow Chart using novaLUM® II-X System


Refer to Current Charm® Operator's Manual for Complete Procedure




Select PLANS button.

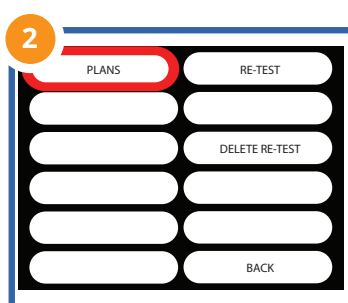

Select PLANS.

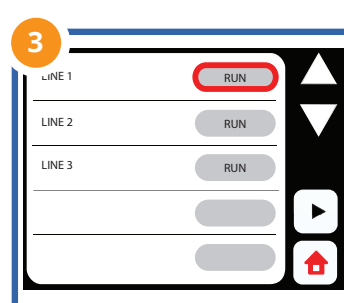

Select RUN next to WaterGiene plan.

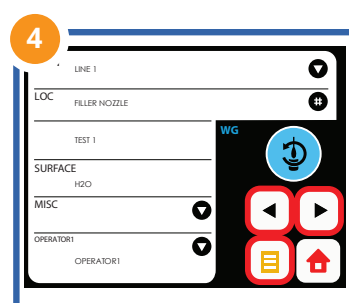

Use Left (◀) or Right (▶) arrows to pick location or use List View (☰) to scroll to find correct location.
Optional: enter operator ID and/or MISC note



Withdraw the swab by gently pulling and twisting the swab handle.



Saturate swab by dipping into liquid for 5 seconds.




Reinsert the swab and twist the handle down. Gently shake side to side 3 times.




Unscrew swab handle completely and shake WaterGiene test 3 times while holding handle so that liquid is forced down into the microtube. Twist handle down again.




Remove black microcap and verify the white tablet is broken up before inserting into the novaLUM II-X system.

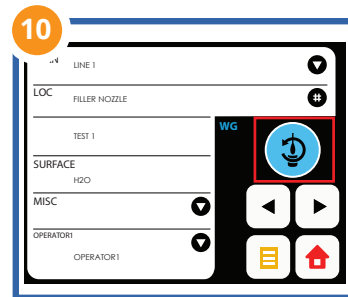

Insert the WaterGiene swab fully into the chamber of the novaLUM II-X system in an upright position. Press (⏻) to start count.

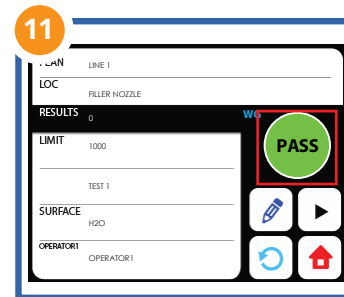

Results in 20 seconds. The result of "Pass" or "Fail" is based on established limit.

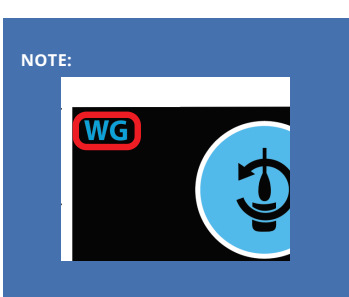

NOTE:
For WaterGiene, the location screen must list WG to the left of the Run Test icon (⏻).