




FieldSwab® Procedural Flow Chart using novaLUM® II-X System


Refer to Current Charm® Operator's Manual for Complete Procedure




Select PLANS button.

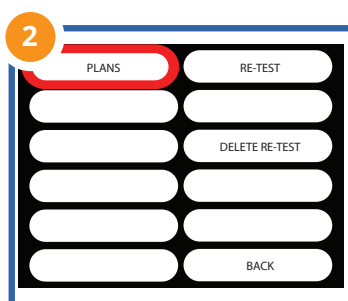

Select PLANS.

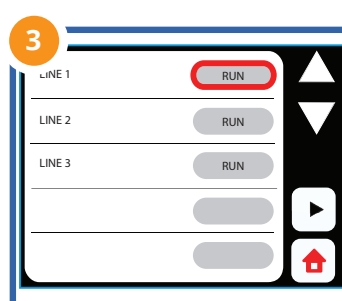

Select RUN.

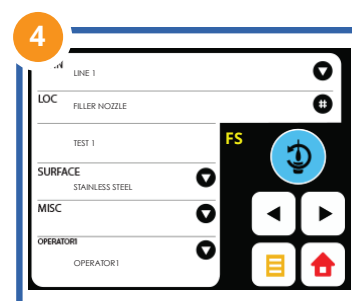

Use Left (◀) or Right (▶) arrows to pick location or use List View (≡) to scroll to find correct location.



Withdraw the swab by gently pulling and twisting the swab handle.




Swab an area 100 cm² (16 in²). Rotate the handle to achieve full swab contact with surface.




Reinsert the swab and twist the handle down. Gently shake side to side 3 times. White tablet will break apart.



Remove black microcap before inserting into the novaLUM II-X system.

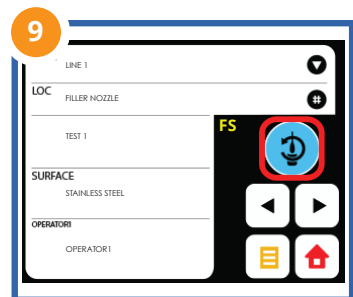

Insert the FieldSwab ATP swab fully into the chamber of the novaLUM II-X system in an upright position. Press (⏻) to start count.

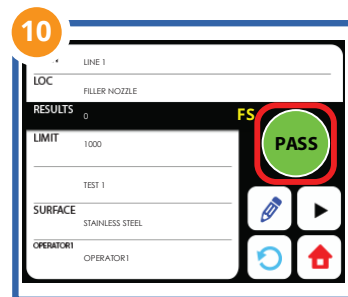

Results in 5 seconds. The result of "Pass" or "Fail" is based on established limit.