




AllerGiene® Procedural Flow Chart using novaLUM® II-X System

Refer to Current Charm® Operator's Manual for Complete Procedure


DO NOT REMOVE BLACK COVER UNTIL AFTER ACTIVATION AND IMMEDIATELY BEFORE INSERTING (ALLERGIENE TEST) INTO THE NOVALUM II-X SYSTEM.




Select PLANS button.

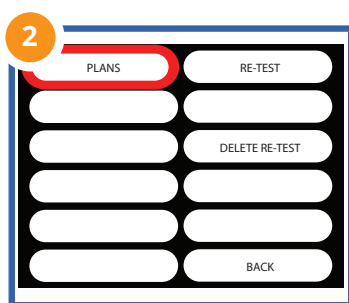

Select PLANS.

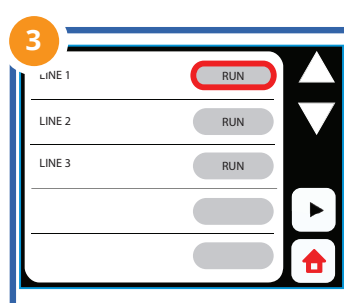

Select RUN.

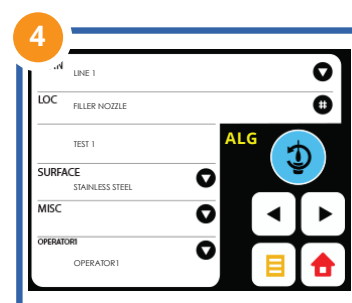

Use Left (◀) or Right (▶) arrows to pick location or use List View (≡) to scroll to find correct location.
Optional: enter operator ID and/or MISC note.




Withdraw the swab by gently pulling and twisting the swab handle.



Swab an area 100 cm² (16 in²). Rotate the handle to achieve full swab contact with surface.




Reinsert the swab and twist the handle down. Gently shake side to side 3 times.

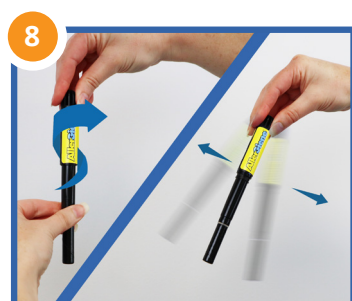


Unscrew swab handle completely and shake AllerGiene 3 times while holding handle so that liquid is forced down into microtube. Twist handle down again.




Remove black microcap and verify the white tablet is broken up before inserting into the novaLUM II-X system.

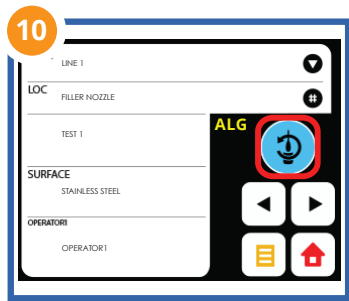

Insert the AllerGiene swab fully into the chamber of the novaLUM II-X system in an upright position. Press (⏻) to start count.

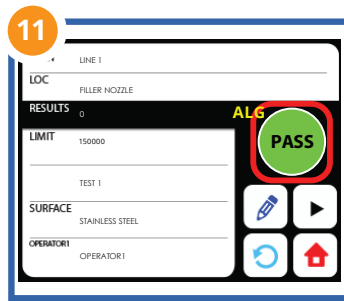

Results in 20 seconds. The result of "Pass" or "Fail" is based on established limit.

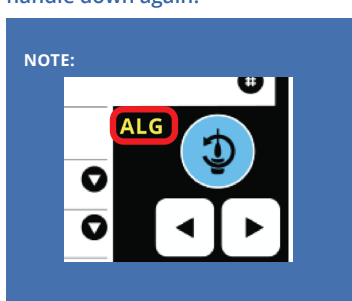

For AllerGiene, the location screen must list ALG to the left of the Run Test icon (⏻).

www.charm.com

T +1.978.687.9200 | F +1.978.687.9216 | E info@charm.com | 659 Andover Street, Lawrence, MA 01843-1032, USA

© 2023 Charm Sciences, Inc. AllerGiene, Charm, and novaLUM are registered trademarks of Charm Sciences, Inc. | MRK-1485-001 Jan-2023

